
When design and function come together...you get FUSION!

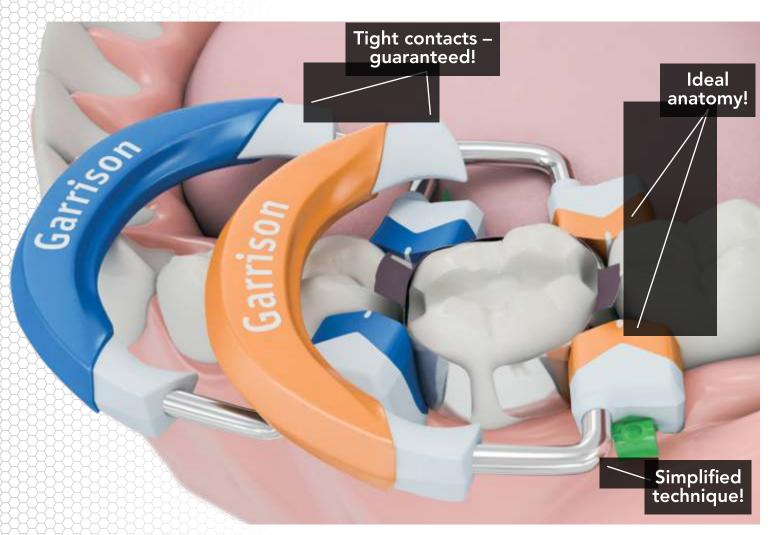

Tight contacts – guaranteed!Superior tooth separation force through advanced engineering and material selection guarantees perfect

Ideal anatomy!

contacts every time.

Soft-Face™ silicone ring tips with enhanced marginal ridge anatomy hug the pre-contoured dead-soft matrices to the tooth for flash free, perfectly contoured restorations.

Simplified technique!

Composi-Tight® 3D Fusion™ is packed with features like easy ring stacking, over-the-wedge placement and Ultra-Grip™ anti-ring slip extensions that make system use intuitive and frustration free.



Only a sectional matrix can deliver results like this!


Tofflemire-style System

Fails to restore proximal anatomy

Thin contact at the marginal ridge

\(\sum_{\text{Large food trap below}} \)

Increased likelihood of fracture, occlusal interference, recurrent caries and periodontal disease

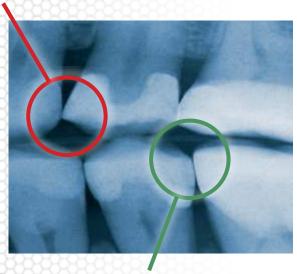

Composi-Tight® Sectional Matrix System



Naturally contoured bands

✓ Anatomically correct contacts

Contacts at the height of contour

Rings produce optimal tooth separation for excellent, tight contacts

Only Composi-Tight® 3D Fusion™ can make the technique this intuitive and frustration free!

1 Place the matrix band:



Dead-soft three dimensionally contoured matrices have placement tabs for easier occlusal placement and Garrison's revolutionary non-stick coating for hasslefree removal.

2 Wedge firmly:



Composi-Tight® 3D Fusion™ Ultra-Adaptive wedges glide in but won't back out. Their soft, flexible exterior readily adapts to interproximal irregularities for a superior gingival seal.

3 Place the ring:



All three of Composi-Tight® 3D Fusion's separator rings will fit right over the wedge for easier placement. With Ultra-Grip™ retention extensions, they will stay firmly in place from the distal of the canine to the widest of molar preparations.



An all-new system designed to handle even the toughest cases

More than 20 years of sectional matrix leadership and an incredible amount of feedback from users all over the globe were poured into the design. All this hard work has resulted in separator rings with significantly improved performance and ease of use.



Stack it up – Tru-Stack[™] technology allows placement of the orange ring (tall) over the blue ring (short) for easier MOD and multiple tooth restorations

Get a grip! All Composi-Tight®
3D Fusion™ separator rings
incorporate Ultra-Grip™
retention extensions that
anchor them securely to
virtually any tooth. These rings
stay where you place them!





Hugs the curves –
Soft-Face™ silicone tips
with enhanced marginal
ridge anatomy hug the
matrix band to the tooth
for flash free restorations

Wider indications, wider smiles!



Dentistry's first Wide Preparation separator ring (green) makes Composi-Tight® 3D Fusion™ the most versatile and user-friendly sectional matrix system available. The Wide Prep ring greatly simplifies what until now has been one of the most challenging posterior composite restorations. Now you can restore these big preps with confidence knowing that you will achieve excellent contact and contour without having to take extra steps.

Matrix band perfection

Composi-Tight® 3D Fusion™ Full Curve Non-Stick Matrices are three dimensionally contoured, ultra-thin and dead-soft to help you recreate ideal interproximal anatomy.

Easy in – easy out! A Grab-Tab[™] centered on the occlusal edge of the band simplifies placement while Slick Bands[™] non-stick coating reduces composite adhesion by an impressive 92%!





Go deep or wide – Improved subgingival extensions and added wraparound length make for the perfect matrix band to use with the green Wide Preparation ring. Now you're ready to tackle even the widest preps!

Composi-Tight®
3D Fusion™
forceps for
precise
control

Composi-Tight® 3D Fusion™ includes dentistry's premier drop-forged ring placement forceps. Regardless of your hand size or strength, Composi-Tight® 3D Fusion™ forceps will give you incredible control for the precise placement of superstrong separator rings over even the widest molars.

Ergonomic for the masses — Whether your hands are big or small, Composi-Tight® 3D Fusion™ forceps give you the leverage you need for precise ring placement.

Garrison



A real softie Composi-Tight® 3D Fusion's Soft-Face™ overmold allows the wedge to do what no other wedge can truly do – actually adapt to interproximal irregularities.

Tough on the inside
The firm inner skeleton allows for
easy insertion and tooth separation
just like a traditional wedge.

Long-term investment – Drop-forged stainless steel prevents leverage-robbing flexion and makes these forceps incredibly durable.

Functionality for all – As an added bonus, Composi-Tight® 3D Fusion™ forceps will work with virtually any ring from any manufacturer.



Only Garrison® can bring you a comprehensive kit like this!

	FX- KFF- 40		FX- KFF - 00		FX- KFF - 14		FX- KFF - 11	
Fusion Rings	Short Tall Wide	01 01 01	Short Tall Wide	01 01 01	Short Tall Wide	01 01 01	Short Tall Wide	01 01 01
Full Curve Non-stick Matrices	4.4mm Gray 6.0mm Red 5.6mm Purple 6.6mm Green 8.7mm Blue TOTAL	10 05 10 10 05 40	4.4mm Gray 6.0mm Red 5.6mm Purple 6.6mm Green 8.7mm Blue TOTAL	20 05 20 20 05 70	4.4mm Gray 6.0mm Red 5.6mm Purple 6.6mm Green 8.7mm Blue TOTAL	20 05 20 20 05 70	4.4mm Gray 6.0mm Red 5.6mm Purple 6.6mm Green 8.7mm Blue TOTAL	20 05 20 20 05 70
Ultra Adaptive Wedges	XS Yellow S Blue M Orange L Green TOTAL	10 10 10 10 40	XS Yellow S Blue M Orange L Green TOTAL	20 20 20 20 20 80	XS Yellow S Blue M Orange L Green TOTAL	20 20 20 20 80	XS Yellow S Blue M Orange L Green TOTAL	20 20 20 20 20 80
Fender Wedges					Assorted	20	Assorted	20
3D Fusion Forceps	1 pc		1 pc		1 pc		1 pc	
3D Fusion Band forceps					1 pc		1 pc	
Multifunction Composite Instrument							1 pc	















