
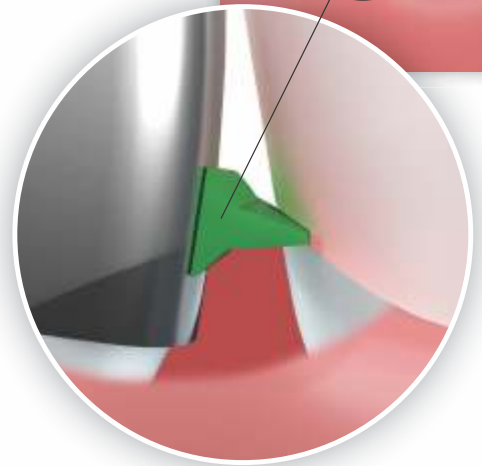


Fusion™ Anterior Matrix System

The all-new **Fusion™ Anterior Matrix System** from Garrison® gives you the materials and techniques to restore many types of anterior lesions with confidence.

For years, our customers have shown us how to use posterior metal matrices to restore anterior teeth. When placed vertically, what was the buccal-lingual curvature of the band creates ideal gingival-incisal curvature for an anterior tooth. We worked with these customers to design a metal matrix band that takes this technique to a whole new level of refinement

Add to that the intelligent design of the new **Fusion™ Anterior Wedges** and you get a complete system. These radically curved wedges hug the bands to the tooth all the way from the facial aspect around to the lingual, helping to free up your hands. Composite placement and sculpting are simplified, allowing you to easily express your artistic nature restoring beautiful smiles for your patients.

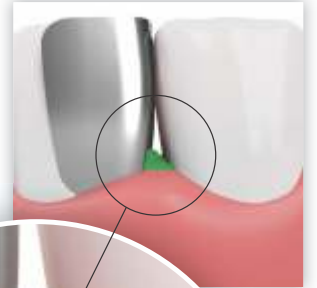

Try the
all-new **Fusion™**
Anterior
Matrix System
today!

Garrison.
Dental Solutions

**The
Social
Six
never
had
it so
good!**

Fusion™ Anterior Matrix System

The **Fusion™ Anterior Matrix** contains a myriad of refinements to help you recreate ideal interproximal anatomy quickly and confidently.



- Only 0.0015" - significantly thinner than strips to facilitate better interproximal contact.
- Firm stainless steel resists deformation both during placement and while compressing composite into the preparation.
- Easy to slide into the sulcus for deeper restorations.
- Ideal anatomic curvature in the gingival-incisal direction AND facial-lingual.
- Far superior to strips when restoring large or deep carious lesions.

Fusion™ Anterior Bands

Refills	Item	Standard		Item	Bulk	
		Qty.	Price		Qty.	Price
Short	AN100-M	50	\$78	AN100	100	\$137
Tall	AN200-M	50	\$78	AN200	100	\$137

The clever features of the **Fusion™ Anterior Wedge** make it a 'must have' for all anterior restorations.



- Radical curvature ensures a firm seal at the cervical margin from facial to lingual.
- Helps maintain ideal anatomy by hugging the band to the tooth.
- Unique 'T' shape sits deeper in the interproximal space to minimize the dreaded 'black triangle.'
- Tapered separator strut allows you to adjust the amount of pressure exerted to control the interproximal contact.
- The cervical interface is half above and half below the sulcus for maximum adaptation with deeper restorations.

Be sure to check out the website at www.garrisonsdental.com

Garrison.
Dental Solutions



1102 Embassy Chambers, 3rd road, Khar West, Mumbai 400052
 +91 22 6127 7777
 sales@nkpatelandsons.com
 +91 99300 06127
 @nkpatelandsons
 /nkpatelandsons